


SPACETOUCH

COSMO PRO

WHOLE BODY
SPECTROMETER

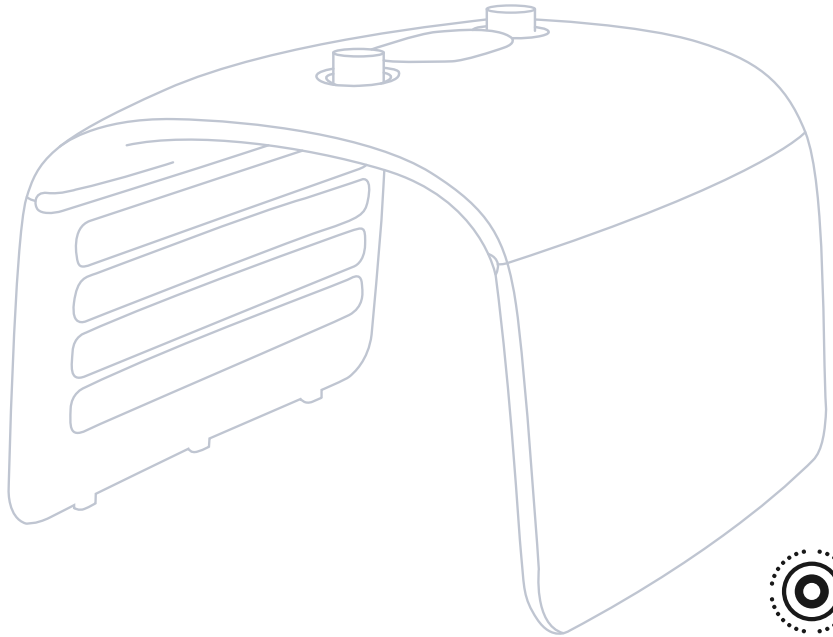

INSTRUCTIONS



COSMO PRO USER GUIDE

HELLO!

Thank you for purchasing Cosmo Pro and welcome to the Spacetouch family!

There is no better investment than investing in yourself! You are now part of a growing community of millions of people who use advanced LED skincare therapy to achieve flawless skin.

As you unbox your new device, take a few moments to review the manual and become familiar with its features before using your Cosmo device. The instructions provide important details on how to safely operate your device.

Spacetouch is committed to providing you with the best skincare experience possible. If you have any questions or concerns, our friendly and knowledgeable customer service team is here to help. You can reach us at Service@Spacetouch.com.

TABLE OF CONTENTS

1	Package Contents	03
2	Cosmo Pro Design	03
3	Device Setup	04
4	Operating Instructions	04
5	Button Overview	04
6	Safety Precautions	09
7	Electrical Precautions	10
8	Health Precautions	11
9	Product Specification	12
10	Troubleshooting	12
11	Warranty Information	13

For the correct operation of this product, please thoroughly read the manual. Retain it for future reference after reading

Important: This manual contains crucial information about device use, safety guidelines, storage recommendations, and warranty details. Ensure you fully understand and follow these safety precautions before operating the device. Failure to do so can result in personal injury or damage, for which we cannot be held responsible.

PACKAGE CONTENTS

Cosmo Pro body device

Adapter

Protective Glasses

ADJUSTABLE DESIGN

We stacked the Cosmo Pro with the latest skin treatment technologies, and made it easy to use at home. With 5 wavelengths, patented serum nano mister, and over 1800 light points, this device illuminates you from head to toe. Its foldable design lets you adjust it to target all body parts.

You can use it in various positions:

- Lying flat
- Lying on your side
- Sitting
- Lying face down



DEVICE SETUP

1. Carefully open the packaging and gently remove the device. Unfold it. Then, place it on a flat, stable surface.
2. Locate the power adapter. Connect its DC end to the device's DC input, ensuring a secure fit.
3. Plug the AC end of the power adapter into a wall outlet. The device will make a "ding" sound upon connecting, indicating it's powered on and in standby mode.
4. The default timer is set to 30 minutes.


OPERATING INSTRUCTIONS

1. Gently tap (short-press) the Power button to turn on the device. You'll hear a "ding" sound, indicating it's been activated successfully.
2. Please avoid pressing or tapping the buttons forcefully.
3. Using the device for sessions lasting between 15 to 30 minutes is recommended.
4. Adjust the session duration as per your requirements.
5. For optimal results, use the device daily.
6. Do not use the device for more than 60 minutes in a single session.

BUTTON OVERVIEW

TOUCH CONTROL SYSTEM

All buttons function through touch control. A "Ding" sound indicates a successful button press. The status of the device is shown using LED indicators.

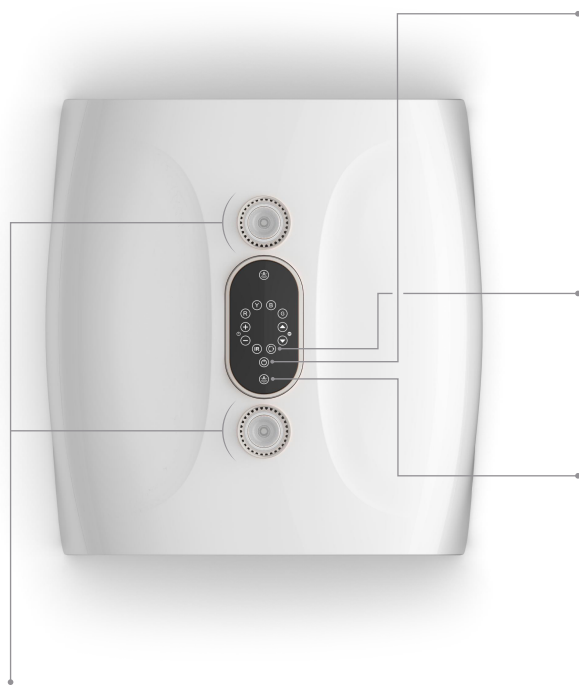

POWER BUTTON (ON/OFF)

- When powered on, the "standby mode" LED indicator will illuminate.
- Tap this button (short press) to turn the device on for a preset duration of 30 minutes.
- Hold the button (long press) to turn off the device and return to "standby mode."

AUTOMATIC COLOR CYCLE BUTTON

- Press once: The buzzer will sound, and the device will begin cycling through red, blue, green, and yellow colors.
- Press again: The color cycling function will be turned off.

NANO AMPOULE MIST BUTTON

- Press this button to activate the mist spray function.
- You'll hear a beep, indicating it's ready to spray.
- You can use the spray function by itself or combine it with other modes.
- Note: This function does not have any adjustable settings.

Important: Before using the mist feature, ensure the two copper points on the bottom of the ampoule connect securely with the main device. This will ensure the mist works correctly.

Cosmo Pro harnesses the power of ultrasonic waves to transform serums and skin care products into a fine mist of nano-sized droplets. These minuscule particles are designed to deeply hydrate your skin and enhance the absorption of the skincare product you're using.

The Cosmo Pro features two mist dispensers: one at the front and one at the back of the device. You have the option to activate either dispenser independently or both simultaneously. You can use the same solution in both or different formulas as recommended by your dermatologist.

For Optimal Results

- We recommend using ampoules from Spacetouch or other trusted sources.
- Ensure the product has a watery consistency. Avoid using thick solutions.
- Do not use concentrated cosmetics or essential oils.
- For even better hydration and anti-aging benefits, combine the use of Cosmo Pro with a facial mask.



TIMER BUTTON

- Once you've connected the adapter, turn on the system using the power switch.
- Adjust the treatment duration by pressing "+" to increase and "-" to decrease.
- Each increment represents 1 minute.
- The timer is set to 30 minutes by default. It can be set to a maximum of 60 minutes and a minimum of 1 minute.


BRIGHTNESS ADJUSTMENT BUTTON

- Use the "↑" button to increase the light intensity of any color mode and the "↓" button to decrease it.
 - There are 5 brightness levels to choose from:
 - Gear 1:** Lowest brightness
 - Gear 5:** Highest brightness
- If Gear 0 is selected, the system is in standby mode.

INFRARED LIGHT BUTTON

(Wavelength: 780-830nm)

- Press this button to activate the near-infrared light function. You will hear a beep as confirmation.
- To turn off the near-infrared light function, simply press the button again.
- Use the brightness adjustment button to change the light intensity.

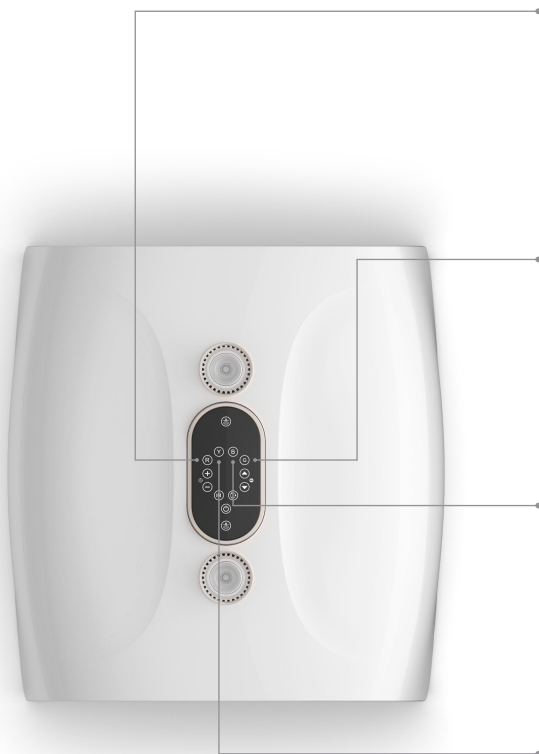

RED LIGHT BUTTON

(Wavelength: 647-760nm)

Press this button to activate the red light function. You will hear a beep as confirmation. To turn off the red light function, simply press the button again.

Use the brightness adjustment button to change the light intensity.

GREEN LIGHT BUTTON

(Wavelength: 492-550nm)

Press this button to activate the green light function. You will hear a beep as confirmation. To turn off the green light function, simply press the button again.

Use the brightness adjustment button to change the light intensity.

BLUE LIGHT BUTTON

(Wavelength: 430-455nm)

Press this button to activate the blue light function. You will hear a beep as confirmation. To turn off the blue light function, simply press the button again.

Use the brightness adjustment button to change the light intensity.

YELLOW LIGHT BUTTON

(Wavelength: 560-620nm)

Press this button to activate the yellow light function. You will hear a beep as confirmation.

To turn off the yellow light function, simply press the button again.

Use the brightness adjustment button to change the light intensity.

SAFETY PRECAUTIONS

GENERAL PRECAUTIONS

- Before using the device, inspect its body and the LED module for any visible damages or defects.
- Before using, ensure there's enough space for the device. Position it away from heaters and other medical equipment.
- Ensure the device is placed in a well-ventilated area. Avoid locations with high temperatures, excessive humidity, or direct sunlight. Make sure the area has the recommended temperature and humidity ranges:
 - Operating Temperature: 0°C to 26°C
 - Storage Temperature: 0°C to 35°C
 - Relative Humidity: Below 80%
- Position the device on a flat surface, ensuring it's level.
- Be cautious when handling the instrument; even slight impacts can cause significant damage. Do not subject it to external shocks, and always handle it carefully.
- Only use official accessories with the machine. Using unapproved accessories or decorations might cause the device to malfunction. Any resulting damages will not be covered by warranty or repaired free of charge.
- Do not attempt to disassemble, repair, or modify the device without proper authorization. Any issues arising from unauthorized disassembly or repair will void the warranty.
- Handle with Care:
 - Avoid dropping, hitting, or placing objects (like mobile phones or magnetic cards) on the machine.
 - Never insert conductive objects, such as screws or wires, into the machine to prevent risks like burns or electric shocks.
 - Minors and those unable to operate or perceive things normally should use the machine under supervision or as advised by a doctor.

- **Cleaning:**
 - Important: Always disconnect the power supply and unplug the power cord before performing any maintenance or cleaning.
 - Turning the power on during these activities can result in equipment malfunctions.
 - Gently wipe the machine with a soft cloth.
 - Avoid using strong chemical agents.
 - Clean tight spots or crevices with a soft brush.
- **Storage:**
 - Ensure that the power is turned off and the power cord is unplugged.
 - Store in a cool, dry, and well-ventilated location.
 - Keep out of reach of children.
 - Avoid direct sunlight, excessive moisture, extreme temperatures, and areas with strong acids, alkalis, oxidants, or heavy dust.

ELECTRICAL PRECAUTIONS

- **Avoid Using Other Heating Devices:** Do not operate this device simultaneously with other heat transfer products. This can lead to malfunction or damage.
- **Safe Distance:** Maintain a safe distance from other electrical appliances. Keep the machine away from strong magnetic sources and flammable items (e.g., gasoline, fireworks).
- **Water and Moisture:** Never pour water on the machine or operate it with wet hands. Ensure that no water enters the device. Moisture can create a high-voltage current that might lead to fires or electric shocks.
- **Power Plug and Sockets:** Refrain from using damaged power plugs or loose sockets. They can pose serious risks.
- **Unplug When Not in Use:** If you're not using the device for an extended period, unplug it. After cleaning, make sure to store it properly. Leaving it plugged in unnecessarily can cause electric shocks, electric leakage, or fires.

- **Emergencies:** If the device ever falls into water, makes an unusual noise, emits a burnt smell, or experiences any other abnormalities, immediately unplug it from the power source.
- **Safety First:** Always prioritize safety and regularly inspect your device for any signs of wear or damage.

HEALTH PRECAUTIONS

- Recommended usage is 15 to 30 minutes per session. This can be adjusted based on individual needs.
- Use 1 to 2 times daily. Do not exceed 60 minutes per session.
- Refrain from using if:
 - You have heart disease, mental illness, or hypertension.
 - You have an infectious skin disease.
 - You suffer from skin inflammation or similar conditions.
 - Even if your condition isn't listed, do not use the device if advised against it by your doctor.
- Never look directly at any light emitted by the machine. When using on your face, always wear protective sunglasses or light-blocking goggles.
- If you have sensitive skin, first test the machine on a small, hidden area of your skin. Only use on larger areas if no allergic reaction occurs.
- Stop using immediately if you experience discomfort or feel unwell after using the product. Before using it again, start with shorter sessions more frequently. Gradually increase the time as your body becomes accustomed.

PRODUCT SPECIFICATION

Product name	Cosmo Pro: Whole body spectrometer
Place of Origin	China
Power Supply	DC24V 12.5A 50HZ/60HZ
LED	Red/Blue/Green/Yellow
Product size	710×690×500 (mm)
Product weight	15kg ±10%

BEFORE REPORTING AN ISSUE

- Check for software updates.
- Determine the potential cause.
- Verify all relevant components are connected.

POWER SUPPLY ISSUES

Loose Adapter Connection: Ensure the adapter is securely connected. Check the power ON/OFF status.

DEVICE STARTUP ISSUES

LED Light Malfunctions: If the LED light fails during a power outage, replace the LED unit. If unsure, contact your supplier or our service center for a replacement.

ACTUATOR PROBLEMS

If you encounter issues with the actuator, please get in touch with your supplier or our service center for assistance.

WARRANTY INFORMATION

Our products come with a lifetime manufacturer warranty, covering any defects under regular use.

Note: This warranty only covers the main device body. Warranty does not cover any misuse of the misting device.

EXCLUSIONS FROM THE FREE WARRANTY

Your device may not be eligible for free service under the following conditions:

- If the device gets damaged due to accidents, mishandling, water damage, or natural disasters (such as fire, salt corrosion, floods, falls, earthquakes, lightning, etc.)
- If the device is modified, repurposed, or disassembled by the user or unauthorized parties.

Remember always to follow the device startup procedures and adhere to the safety guidelines in this manual to maintain your warranty coverage.

IMPORTANT NOTES

- We take pride in our strict quality control processes. To avail of free maintenance, please provide the following warranty certificate.
- During the warranty period, we offer free maintenance for any malfunctions or faults that arise during standard use.
- Suppose the equipment fails due to reasons such as using a non-original power supply, misuse, unauthorized disassembly, or any man-made damage. In that case, we reserve the right to deny free maintenance and may charge for parts and services.
- A second copy of this warranty certificate will not be issued. Please ensure it's kept safe and secure.

WARRANTY CERTIFICATE

Model _____

Product Name _____

Date of Purchase _____

Corporate Name _____

Mailing Address _____

Full Name _____

Please fill in the above details and keep this certificate safe
for any future warranty claims.

COSMO PRO

WHOLE BODY SPECTROMETER

Designed in California
By Spacetouch
Manufactured by Pandmedic inc
4290 Cameron St unit 7
Las Vegas NV B9103

Customer Support
Phone: #1 833-801-0251
Email: Service@spacetouch.com

