SPACETOUCH

ATLAS

THE SMART WELLNESS DEVICE





ATLAS USER GUIDE

HELLO!

Thank you for purchasing Atlas and welcome to the Spacetouch family!

There is no better investment than investing in yourself! You are now part of a growing community of millions of people who use advanced LED skincare therapy to achieve flawless skin.

As you unbox your new device, take a few moments to review the manual and become familiar with its features before using your Cosmo device. The instructions provide important details on how to safely operate your device.

Spacetouch is committed to providing you with the best skincare experience possible. If you have any questions or concerns, our friendly and knowledgeable customer service team is here to help. You can reach us at Service@Spacetouch.com.

TABLE OF CONTENTS

1	Overview	03
2	Customizable Treatments	03
3	Features	04
4	Voice Commands Guide	05
5	Touchscreen Controller	08
6	Pre-Programmed Treatments	09
7	Mannual Mode	32
8	Entertainment	34
9	Precautions	35
10	Product Specification	38
11	Troubleshooting	38
12	Warranty Information	39

We have made it even easier to use our best-selling full-body wellness device for both spa and home use.

The **Atlas** is equipped with **11 preset advanced programs** that can be activated with a touch of a button or by voice command. These programs cater to various needs, including **pain relief**, **detox**, **muscle recovery**, and **stress relief**. The Atlas is large enough to fit any area of the body yet light enough to carry.

Its unique design allows for easy folding and storage in a protective case, making it a hit among both home users and spa owners who wish to add advanced skincare technologies but are limited in space.

CUSTOMIZABLE TREATMENTS

In addition to the preset programs, you have access to a **Manual Mode** for those seeking even more customized treatments. The Manual Mode allows you to control:

- Light Spectrum
- Light Intensity
- Heat Level
- Treatment Duration

PATENTED SERUM MISTER

Atlas is equipped with two patented serum misters, which automatically distribute skincare products to the treated area. The combination of nano mist with light therapy has shown an increase of 63% in results compared to separate treatments.

EFFORTLESS VOICE CONTROL

For even easier control, you can manage your treatment by simply talking to Atlas without the need to stop the session. Using voice commands, you can control programs, adjust heat, add or reduce time, start or stop the mist, and start or stop programs. This exciting feature provides a seamless wellness experience for your treatment.

SEAMLESS WELLNESS EXPERIENCE

Wellness is all about convenience and feeling better from the inside out. While treating your skin, enjoy a seamless experience with Atlas by connecting via Bluetooth to play your favorite music or podcasts during your treatment.

KEY FEATURES

11 Wellness Programs

- Suitable for Home and Spa Use
- Treats Any Body Area

Noice Control

Bluetooth Connection

Built-in Speakers/Sound System

Touchscreen Interface

FDA Made in an FDA-Approved Facility

VOICE COMMANDS GUIDE

ACTIVATION

• "Hey Atlas" - Activates the device.

LIGHT CONTROL

- "Lights On" Turns on the light.
- "Lights Off" Turns off the light.

PROGRAM SELECTION

- "Mode One" Activates Body and Pain Relief program.
- "Mode Two" Activates Anti-Aging & Firming program.
- "Mode Three" Activates Detox & Glow program.
- "Mode Four" Activates Skin Discoloration program.
- "Mode Five" Activates Cellulite & Stretch Marks program.
- "Mode Six" Activates Muscle Recovery program.
- "Mode Seven" Activates Stress Relief program.
- "Mode Eight" Activates Dry Skin program.
- "Mode Nine" Activates Acne & Oil Control program.
- "Mode Ten" Activates Rejuvenation & Repair program.
- "Mode Eleven" Activates Body Slimming & Contouring program.

LIGHT SPECTRUM CONTROL

- "Red Light" Activates red light.
- "Blue Light" Activates blue light.
- "Yellow Light" Activates yellow light.
- "Green Light" Activates green light.
- "Infrared Light On" Turns on infrared light (can be used with other lights).
- "Infrared Light Off" Turns off infrared light.

NANO SERUM MISTER CONTROL

- "Mist On" Turns on the nano serum mister.
- "Mist Off" Turns off the nano serum mister.

MUSIC CONTROL

- "Volume Up" Increases the volume.
- "Volume Down" Decreases the volume.

ENERGY CONTROL

Adjusts the light and heat strength of the treatment:

- "Level One" Sets treatment strength to level one.
- "Level Two" Sets treatment strength to level two.
- "Level Three" Sets treatment strength to level three.

- "Level Four" Sets treatment strength to level four.
- "Level Five" Sets treatment strength to level five.

TIME CONTROL

Allows you to control the treatment time in 10-minute increments:

- "Ten Minutes" Sets treatment timer for ten minutes.
- "Twenty Minutes" Sets treatment timer for twenty minutes.
- "Thirty Minutes" Sets treatment timer for thirty minutes.
- "Forty Minutes" Sets treatment timer for forty minutes.
- "Fifty Minutes" Sets treatment timer for fifty minutes.
- "Sixty Minutes" Sets treatment timer for sixty minutes.


TOUCHSCREEN CONTROLLER

The Atlas touchscreen provides professional preset treatments designed by skincare experts to help you achieve your desired results with the touch of a button, all in the convenience of your home. It's a great addition to your spa, as you won't need an esthetician to activate the device—your customers can self-service.




PRE-PROGRAMMED TREATMENTS

The Capsule is equipped with a responsive touchscreen interface on its exterior. Simply tap to browse and select from the available program options.

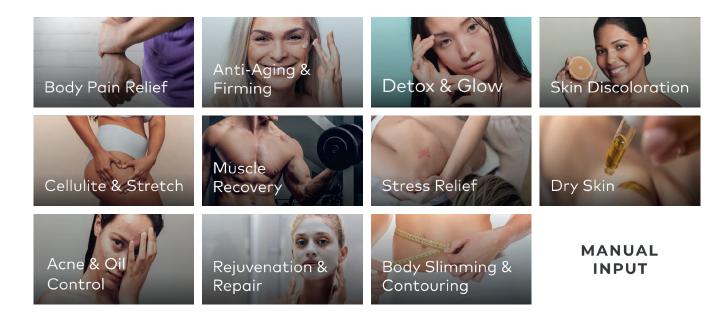


The Body Pain Relief Program is designed to alleviate discomfort, reduce inflammation, and soothe tight muscles through a combination of therapeutic heat, deep-penetrating infrared light, and red LED therapy. This synergistic approach promotes muscle relaxation, enhances circulation, and supports the body's natural healing mechanisms, providing relief from chronic pain, muscle stiffness, and joint soreness.

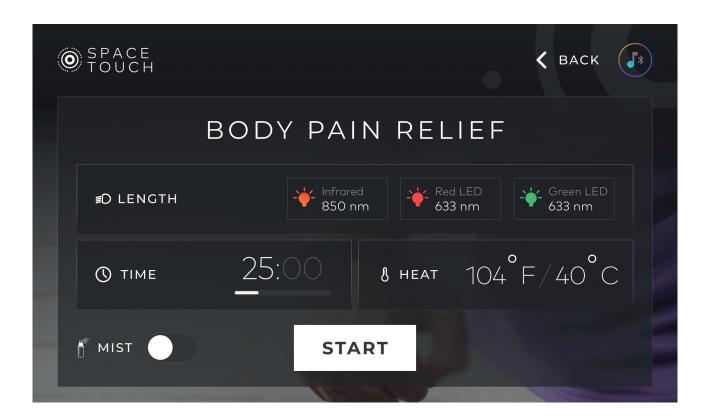





Experience the pinnacle of skin rejuvenation with the Anti-Aging & Firming Program. Combining warmth to enhance product absorption, infrared and red LED light to stimulate collagen production, and yellow light to improve skin tone, this program targets all signs of aging. It visibly reduces wrinkles, improves skin elasticity, and leaves the skin looking firmer and more youthful.

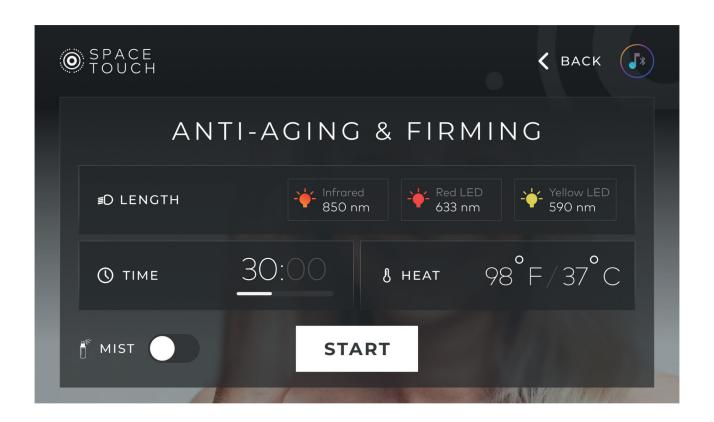




Unlock your skin's natural radiance with our Detox & Glow Program. This luxurious treatment uses heat to detoxify the skin, red LED light to stimulate cellular renewal, and yellow light to brighten the complexion. It's perfect for anyone seeking a noticeable, natural glow and healthier, more vibrant skin.

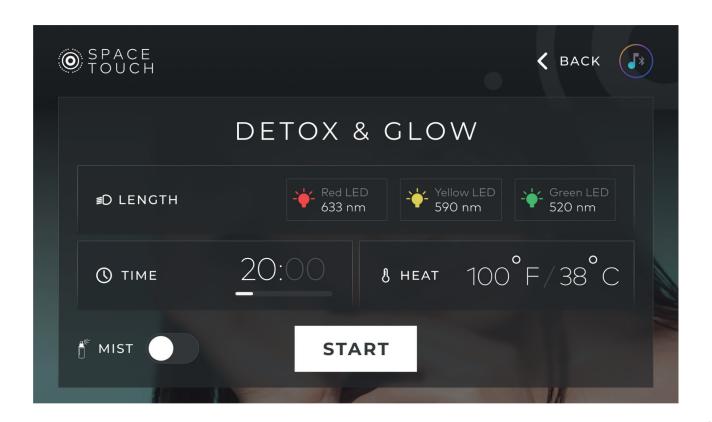

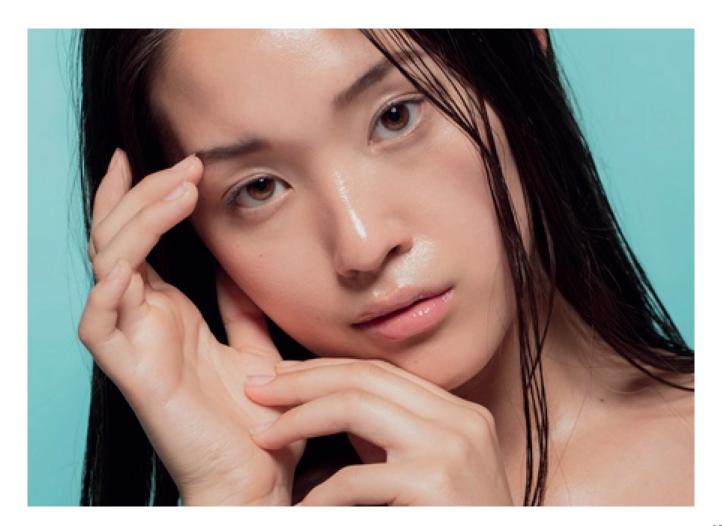



Target skin discoloration with precision using our Skin Discoloration Program. This program combines red and blue LED light therapy to address sun spots, pigmentation issues, and vitiligo by promoting an even skin tone, reducing the appearance of discoloration, and supporting skin health. Achieve clearer, more uniform skin with each session.

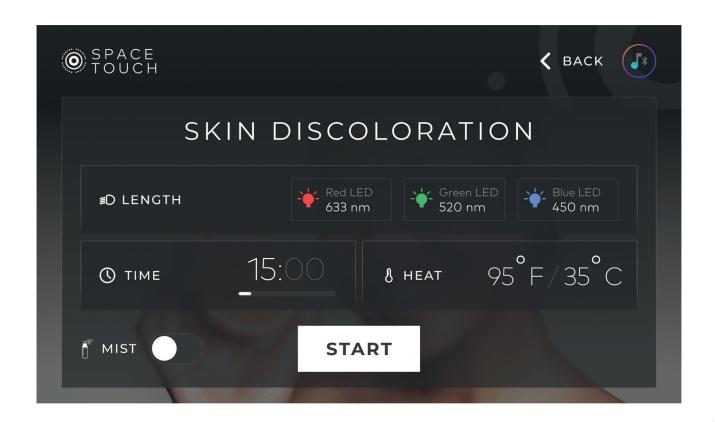




Combat cellulite and stretch marks with our targeted program, designed to smooth skin texture and enhance elasticity. The combination of heat, infrared, and red LED light therapy works to break down fat cells, boost collagen production, and improve the appearance of cellulite and stretch marks for firmer, more toned skin.

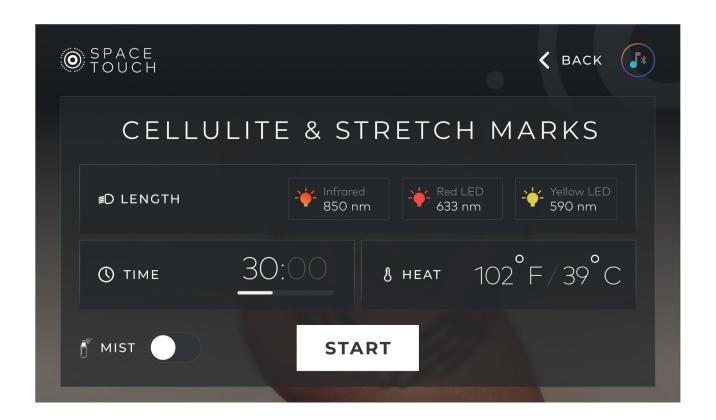

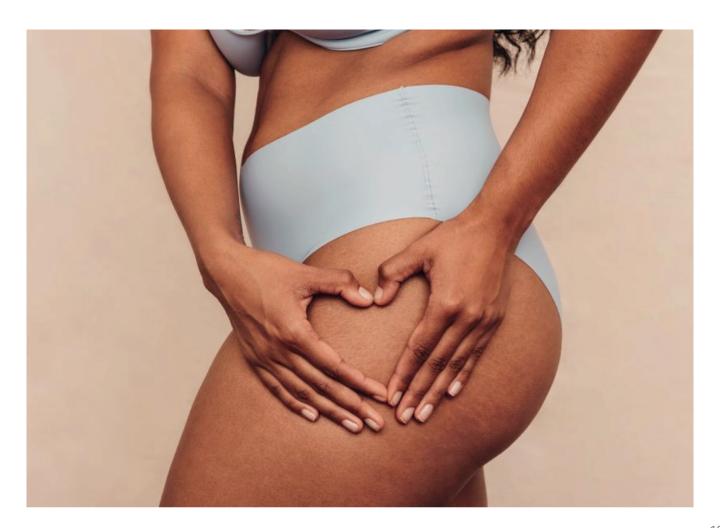


Accelerate muscle recovery and reduce soreness with this powerful program. The deep penetration of infrared light, coupled with therapeutic heat, helps to relax muscles, enhance circulation, and support the body's recovery process after physical activity or strain.

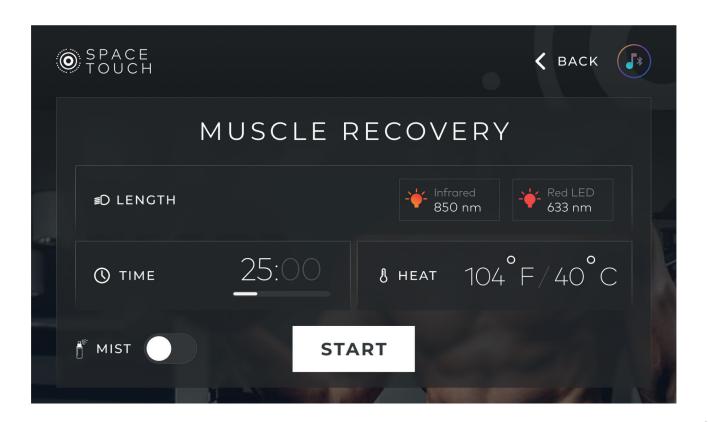




Immerse yourself in a soothing oasis of calm with our Stress Relief Program. Gentle heat and blue LED light therapy, accompanied by relaxing spa music, work together to calm the mind, reduce stress levels, and promote a sense of well-being. Escape the hustle and bustle and emerge refreshed and rejuvenated.

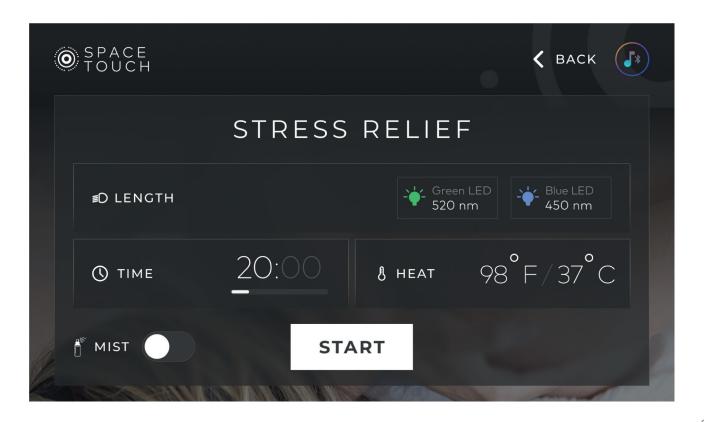




Hydrate and nourish dry skin with our specialized program. The combination of gentle heat and yellow light therapy boosts moisture retention, enhances barrier function, and leaves the skin feeling soft, supple, and rejuvenated. Say goodbye to dryness and welcome a hydrated, glowing complexion.

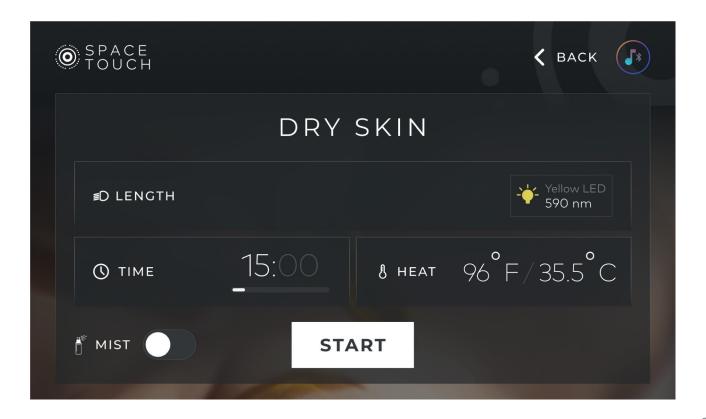





Achieve clear, balanced skin with our Acne & Oil Control Program. Utilizing blue LED light to combat acne-causing bacteria and yellow light to soothe inflammation, this program effectively reduces acne breakouts, controls oil production, and promotes a healthier complexion.

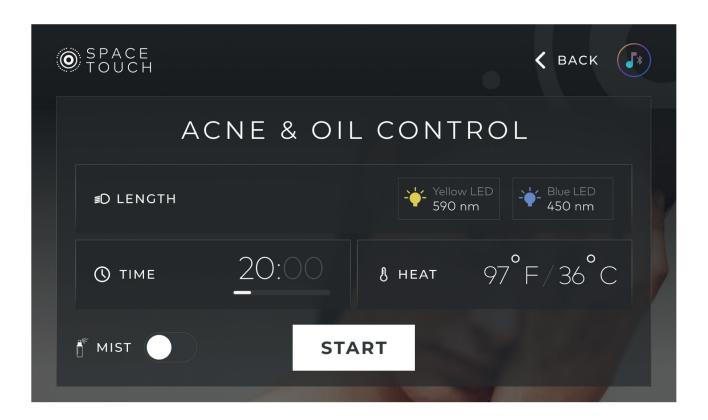




Revitalize your skin with our Rejuvenation & Repair Program, designed to enhance collagen production, accelerate healing, and promote deep skin rejuvenation. The synergy of heat, infrared, and red LED light therapy targets fine lines, wrinkles, and restores a youthful, radiant complexion.

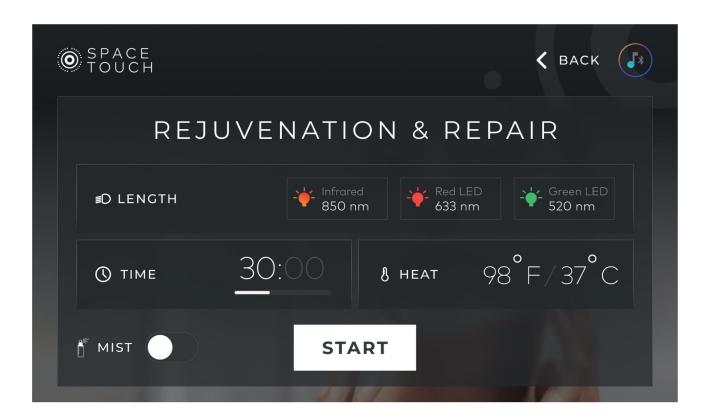

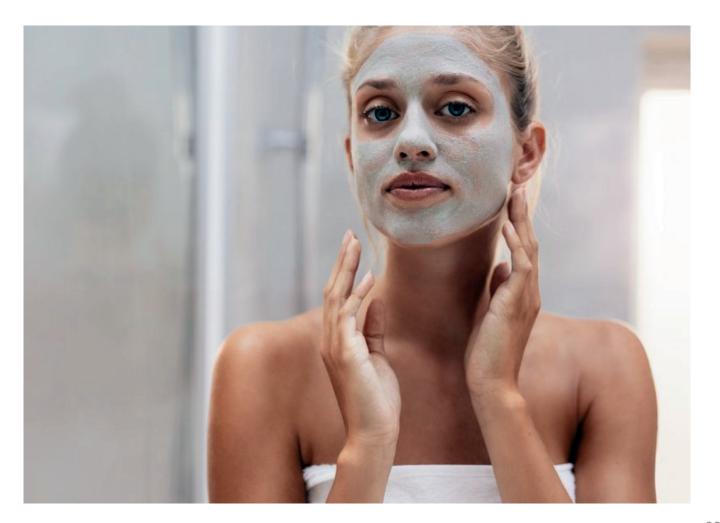

This comprehensive approach warms and prepares the skin, enhances circulation, and promotes detoxification, aiming to refine your silhouette, improve skin tone, and texture. Experience a revitalized, more contoured you with this unique treatment in the Capsule, designed for a firmer, healthier-looking body. Start your transformation today.

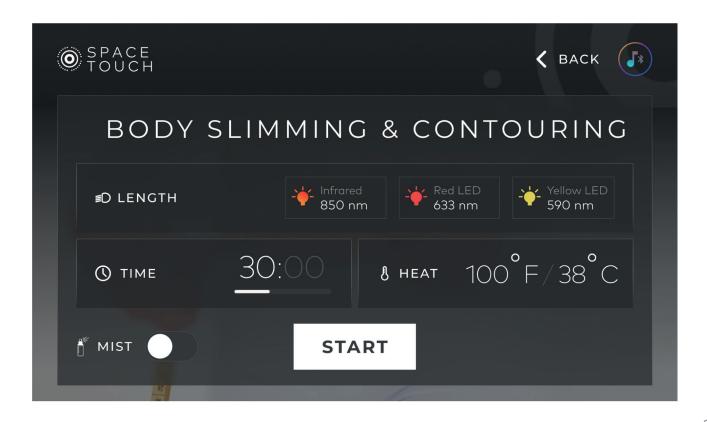



MANNUAL MODE

Think of your body as a rechargeable battery. Just as batteries need power to function, so do we. We constantly interact with our surroundings, drawing energy and vitality.

One of the significant ways our bodies derive this energy is from light. Low-Level Light Therapy (LLLT) harnesses the power of specific light wavelengths, delivered through Light Emitting Diodes (LEDs), to supercharge our body's natural functions.

How Does LLLT Work?

LLLT utilizes specific light wavelengths that our cells readily absorb and convert into energy. As these light beams penetrate our layers of skin, they not only prompt increased circulation but also stimulate the release of Nitric Oxide, our body's natural vasodilator. This energy surge initiates various cellular processes, from healing to regeneration and even the production of essential proteins like collagen.

With its targeted light therapy, the Atlas takes this natural rejuvenation and healing ability of our skin and amplifies it, resulting in enhanced vitality and well-being.





LIGHT WAVELENGTHS AND THEIR BENEFITS

RED LIGHT (633NM)



LED Count: 4,488 LEDs

Depth of Penetration: Up to 10mm

Key Benefits:

- Anti-aging: Boosts collagen production to enhance skin elasticity.
- Improved Blood Circulation: Enhances blood flow to the skin.

When to Use: For addressing signs of aging, improving skin elasticity, and promoting healthy blood circulation.

INFRARED LIGHT (850NM)



LED Count: 2.376 LEDs

Depth of Penetration: Up to 30mm (This offers the deepest penetration among the lights)

Key Benefits:

- Cellulite Treatment: Reduces inflammation to treat underlying fat structures.
- Wound Healing: Accelerates the healing process.
- Rejuvenation: Promotes cellular renewal and health.

When to Use: For cellulite treatment, promoting faster wound recovery and cellular rejuvenation.

BLUE LIGHT (450NM)



LED Count: 2.376 LEDs

Depth of Penetration: Surface level, up to 2mm

Key Benefits:

- Acne Management: Targets and eliminates acne-causing bacteria.
- Oil Control: Helps regulate sebum production.
- Soothing Effect: Calms sensitive and irritated skin.

When to Use: For managing acne, controlling excessive skin oil, and calming skin sensitivity.

YELLOW LIGHT (590NM)



LED Count: 2.376 LEDs

Depth of Penetration: Up to 5mm


Key Benefits:

- Skin Correction: Reduces skin redness and evens out tone.
- Detoxifying Effect: Improves lymphatic flow to detoxify the skin.
- Radiance Boost: Amplifies skin's natural glow and luminosity.

When to Use: For correcting skin discoloration, detoxifying the skin, and boosting overall skin radiance.

ENTERTAINMENT

Enhance your relaxation experience with seamless entertainment control. Connect your phone via Bluetooth to the integrated SpaceTouch music system and enjoy your favorite tunes directly from the device.

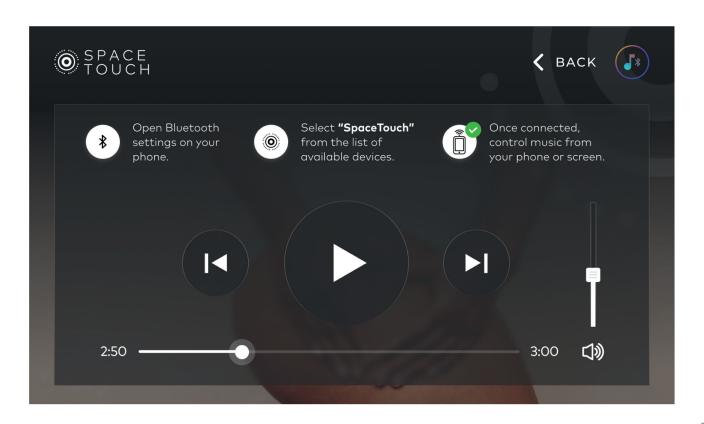

SAFETY PRECAUTIONS

GENERAL PRECAUTIONS

- Before using the device, inspect its body and the LED module for any visible damages or defects.
- Before using, ensure there's enough space for the device. Position it away from heaters and other medical equipment.
- Ensure the device is placed in a well-ventilated area. Avoid locations with high temperatures, excessive humidity, or direct sunlight. Make sure the area has the recommended temperature and humidity ranges:
 - Operating Temperature: 0°C to 26°C
 - Storage Temperature: 0°C to 35°C
 - Relative Humidity: Below 80%
- Position the device on a flat surface, ensuring it's level.
- Be cautious when handling the instrument; even slight impacts can cause significant damage. Do not subject it to external shocks, and always handle it carefully.
- Only use official accessories with the machine. Using unapproved accessories or decorations might cause the device to malfunction. Any resulting damages will not be covered by warranty or repaired free of charge.
- Do not attempt to disassemble, repair, or modify the device without proper authorization. Any issues arising from unauthorized disassembly or repair will void the warranty.
- Handle with Care:
 - Avoid dropping, hitting, or placing objects (like mobile phones or magnetic cards) on the machine.
 - Never insert conductive objects, such as screws or wires, into the machine to prevent risks like burns or electric shocks.
 - Minors and those unable to operate or perceive things normally should use the machine under supervision or as advised by a
 doctor.

Cleaning:

- Important: Always disconnect the power supply and unplug the power cord before performing any maintenance or cleaning.
- Turning the power on during these activities can result in equipment malfunctions.
- Gently wipe the machine with a soft cloth.
- Avoid using strong chemical agents.
- Clean tight spots or crevices with a soft brush.

• Storage:

- Ensure that the power is turned off and the power cord is unplugged.
- Store in a cool, dry, and well-ventilated location.
- Keep out of reach of children.
- Avoid direct sunlight, excessive moisture, extreme temperatures, and areas with strong acids, alkalis, oxidants, or heavy dust.

ELECTRICAL PRECAUTIONS

- **Avoid Using Other Heating Devices:** Do not operate this device simultaneously with other heat transfer products. This can lead to malfunction or damage.
- Safe Distance: Maintain a safe distance from other electrical appliances. Keep the machine away from strong magnetic sources and flammable items (e.g., gasoline, fireworks).
- Water and Moisture: Never pour water on the machine or operate it with wet hands. Ensure that no water enters the device. Moisture can create a high-voltage current that might lead to fires or electric shocks.
- Power Plug and Sockets: Refrain from using damaged power plugs or loose sockets. They can pose serious risks.
- **Unplug When Not in Use:** If you're not using the device for an extended period, unplug it. After cleaning, make sure to store it properly. Leaving it plugged in unnecessarily can cause electric shocks, electric leakage, or fires.

- **Emergencies:** If the device ever falls into water, makes an unusual noise, emits a burnt smell, or experiences any other abnormalities, immediately unplug it from the power source.
- Safety First: Always prioritize safety and regularly inspect your device for any signs of wear or damage.

HEALTH PRECAUTIONS

- Recommended usage is 15 to 30 minutes per session. This can be adjusted based on individual needs.
- Use 1 to 2 times daily. Do not exceed 60 minutes per session.
- Refrain from using if:
 - You have heart disease, mental illness, or hypertension.
 - You have an infectious skin disease.
 - You suffer from skin inflammation or similar conditions.
 - Even if your condition isn't listed, do not use the device if advised against it by your doctor.
- Never look directly at any light emitted by the machine. When using on your face, always wear protective sunglasses or light-blocking goggles.
- If you have sensitive skin, first test the machine on a small, hidden area of your skin. Only use on larger areas if no allergic reaction occurs.
- Stop using immediately if you experience discomfort or feel unwell after using the product. Before using it again, start with shorter sessions more frequently. Gradually increase the time as your body becomes accustomed.

PRODUCT SPECIFICATION

Product name	Atlas: The Smart Wellness Device
Place of Origin	China
Power Supply	DC24V 12.5A 50HZ/60HZ
LED	Red/Blue/Green/Yellow
Product size	710×690×500 (mm)
Product weight	15kg ±10%

BEFORE REPORTING AN ISSUE

- Check for software updates.
- Determine the potential cause.
- Verify all relevant components are connected.

POWER SUPPLY ISSUES

Loose Adapter Connection: Ensure the adapter is securely connected. Check the power ON/OFF status.

DEVICE STARTUP ISSUES

LED Light Malfunctions: If the LED light fails during a power outage, replace the LED unit. If unsure, contact your supplier or our service center for a replacement.

ACTUATOR PROBLEMS

If you encounter issues with the actuator, please get in touch with your supplier or our service center for assistance.

WARRANTY INFORMATION

Our products come with a lifetime manufacturer warranty, covering any defects under regular use.

Note: This warranty only covers the main device body. Warranty does not cover any misuse of the misting device.

EXCLUSIONS FROM THE FREE WARRANTY

Your device may not be eligible for free service under the following conditions:

- If the device gets damaged due to accidents, mishandling, water damage, or natural disasters (such as fire, salt corrosion, floods, falls, earthquakes, lightning, etc.)
- If the device is modified, repurposed, or disassembled by the user or unauthorized parties.

Remember always to follow the device startup procedures and adhere to the safety guidelines in this manual to maintain your warranty coverage.

IMPORTANT NOTES

- We take pride in our strict quality control processes. To avail of free maintenance, please provide the following warranty certificate.
- During the warranty period, we offer free maintenance for any malfunctions or faults that arise during standard use.
- Suppose the equipment fails due to reasons such as using a non-original power supply, misuse, unauthorized disassembly, or any man-made damage. In that case, we reserve the right to deny free maintenance and may charge for parts and services.
- A second copy of this warranty certificate will not be issued. Please ensure it's kept safe and secure.

WARRANTY CERTIFICATE

Model	Product Name
Date of Purchase	Corporate Name
Mailing Address	
Full Namo	

Please fill in the above details and keep this certificate safe for any future warranty claims.

ATLAS

THE SMART WELLNESS DEVICE

Designed in California By Spacetouch Manufactured by Pandmedic inc 5470 Wynn Road, Suite 100, Las Vegas, NV 89118

Customer Support Phone: #1 833-801-0251

Email: Service@spacetouch.com

