

SPACETOUCH

FULL BODY DEVICE



LUNA USER GUIDE

HELLO!

Thank you for purchasing Luna and welcome to the Spacetouch family!

There is no better investment than investing in yourself! You are now part of a growing community of millions of people who use advanced LED skincare therapy to achieve flawless skin.

As you unbox your new device, take a few moments to review the manual and become familiar with its features before using your Luna device. The instructions provide important details on how to safely operate your device.

Spacetouch is committed to providing you with the best skincare experience possible. If you have any questions or concerns, our friendly and knowledgeable customer service team is here to help. You can reach us at Service@Spacetouch.com.

TABLE OF CONTENTS

1	What Is Luna?	04
2	How Does It Work?	05
3	Technical Specifications	08
4	Compliance	08
5	Package Contents	09
6	How to Use Luna	10
7	Treatment Index	15
8	Contraindications, Warnings, and Precautions	32
9	Light Sensitivity Test	33
10	Limited Lifetime warranty	34
11	Register Your device	35



WHAT IS LUNA?

Luna is a light therapy device that combines UVB and LED light in one easy-to-use device. Designed by plastic surgeons and dermatologists, Luna incorporates multiple technologies in one treatment to target a wide array of concerns, including bone density, mood and sleep disorders, jaundice, skin-related cancer symptoms, and skin disorders such as psoriasis or vitiligo. Luna offers 7 modes of treatment- red, blue, green, yellow, orange, violet, and calcium supplementation infrared wavelengths. All treatments can be used on the face, neck, and body. This state-of-the-art device is powered by NASA's revolutionary infrared technology that focuses on the long-term health of your skin. Luna is rechargeable and travel-friendly, meaning that it provides powerful results wherever you are, allowing you to receive your LED light treatments from the comfort of your home.

HOW DOES IT WORK?

Light therapy is a non-invasive treatment that utilizes light at precise wavelengths to penetrate the skin's surface without causing any damage. The light stimulates natural chemical reactions in the skin cells, triggering the rejuvenation and healing processes. Developed by leading plastic surgeons and dermatologists, Luna is a home-use light device that effectively addresses a plethora of various skin concerns via its 7 available modes of treatment:

RED LIGHT SETTING

The red light setting is designed to reduce signs of aging. From the age of 22 onwards, our skin loses between 1-2% of its collagen every year. Collagen plays a critical role as the "glue" between our skin layers, imparting a plump appearance to our skin. As the collagen dwindles, our skin becomes drier and less supple, leading to the emergence of visible lines and wrinkles on the surface of the skin. Furthermore, factors like sleep deprivation, dehydration, sun exposure, and smoking can exacerbate these aging effects. The red light setting is designed to assist the reproduction of cells, collagen, and fiber and reaches the hypodermis level, the deepest and thickest layer of the skin. Consistent use of healthy lifestyle choices and red light LED treatment can effectively halt the reduction of collagen as we age and in some cases, even boost the body's natural collagen production. Therefore, it is crucial to start utilizing red LED light therapy as early as possible.

BLUE LIGHT SETTING

The blue light setting is designed to improve skin texture and reduce sebaceous hyperplasia (enlarged oil glands). It can help with reducing the appearance of actinic keratoses (sun spots), acne, and even acne scarring. This

This treatment can be used in tandem with other products as part of an acne treatment plan. In addition to its cosmetic benefits, blue light therapy has been shown to promote the synthesis of protein and bone collagen, and it can also help with whitening skin and teeth. Blue light therapy is also effective in managing neonatal jaundice, skin conditions such as allergy-based and oily skin, and the relaxation of nerves. It can also provide pain relief and is used for the clinical treatment of goiter, chronic heart disease, insomnia, and menopausal syndrome.

GREEN LIGHT SETTING

The green light setting is designed to stimulate collagen production, reduce inflammation, and improve energy levels. In addition to its cosmetic benefits, green light therapy has been used clinically to treat various health conditions. It is effective in improving symptoms of asthma, relieving climacteric syndrome, and treating nerve numbness and kidney diseases. It has also been used to manage diabetes and digestive tract diseases, making it a valuable treatment option for those seeking to improve their overall health.

YELLOW LIGHT SETTING

The yellow light setting in light therapy offers a unique combination of the benefits of both red and green light. The treatment is designed to reduce pigmentation caused by sun damage, age spots, redness, spider veins, and rosacea. However, yellow light therapy can also accelerate wound healing and improve the healing of ulcers. In addition to its skin-related benefits, yellow light therapy has been shown to enhance the secretion of adrenaline, stimulate the heart clinically, and dilate the trachea during asthma.

ORANGE LIGHT SETTING

The orange light setting in light therapy offers a unique combination of the benefits of both red and yellow light. This makes it an excellent choice for improving the appearance of fine lines and wrinkles, boosting collagen production, reducing inflammation, and stimulating circulation. In addition, orange light has been shown to bring vitality to dull and lifeless complexions. Beyond its cosmetic benefits, orange light therapy can also accelerate the metabolism, increase appetite, and effectively treat emaciation and anemia. It has also been used clinically to treat general heart disease, particularly weak heartbeats.

VIOLET LIGHT SETTING

The violet light setting combines the skin-clearing benefits of blue light with the anti-aging benefits of red light. This makes it an excellent choice for patients looking to address both acne and aging concerns. Violet LED light has been shown to help destroy acne-causing bacteria while also promoting collagen production, resulting in a clearer, younger-looking complexion. In addition to its cosmetic benefits, purple LED therapy is also beneficial for cellular oxygenation and regeneration, which promotes enhanced skin fitness and vitality. This therapy has also been shown to stimulate the secretory system, drain lymph glands, and effectively treat skin and mucosal ulcers.

CYCLE MODE

The cycle mode is designed to provide the ultimate in-skin treatment by combining all LED light modes. This is a fantastic option for those wanting to target multiple concerns in a single treatment and improve general health and well-being. This mode is set automatically to 30 minutes, however, you can manually adjust the timer according to

your needs or your skin specialist's recommendation.

● **CALCIUM SUPPLEMENTATION LIGHT THERAPY**

Luna is the only home-use device to feature calcium supplementation light therapy. Calcium supplementation is excellent for the promotion of good cholesterol, increasing bone mineral density, the treatment of osteoporosis, cartilage diseases, dysplasia, and prevention of complications caused by abnormal metabolism of calcium and phosphorus ions. Additionally, calcium intensely nourishes the skin, softens it effectively, and promotes the strengthening of the skin barrier to help make the skin denser.

TECHNICAL SPECIFICATIONS

Luna weight: 3.0KG / 6.6 LBS

Product measurement: 54 X 53.5 X 29CM / 21.2 X 21.06 X 11.4IN

Input power: 100-240V 50/60 HZ

COMPLIANCE

All Spacetouch LED devices are approved by the U.S. Food and Drug Administration (FDA) and comply with the international standard of basic safety for medical devices IEC 60601 and with the relevant collateral and particular standards for non-laser light source medical equipment for home use IEC 60825-1 and IEC 62471 Photobiological Safety of Lamps and Lamp Systems. Spacetouch devices also conform to the application of risk management to medical devices ISO 14971, and biological evaluation of medical devices in tests for invitro cytotoxicity

ISO 10993-5 and irritation and skin sensitization ISO 10993-10.

PACKAGE CONTENTS

User manual guide

VIP warranty card

Device touch panel

Device side panel (x2)


Power adapter

Protective goggles

HOW TO USE LUNA

TURNING THE DEVICE ON

To start using Luna, you need to assemble the 3 device pieces: the touch panel and 2 side devices. Begin by connecting the first side piece to the touch panel, aligning the 3 tabs of the mounting brackets under the 2 rectangular spaces at the end of the touch panel. Then, slide and press firmly against the panels until the tabs engage. Then repeat the process on the other side of the touch panel with the other side device. Once the device is assembled, plug it into the power source using the cable provided, taking care to only use the original cable as other cables may damage the device. Plug the power source into an electric outlet. Once the power is connected, Luna will beep to indicate that it is ready to use. Please, note that the text on the touch panel will not appear until the device is powered on. To operate, tap the power button and wait for the sound as a response, indicating a successful operation. With these steps completed, you are now ready to select your treatment!

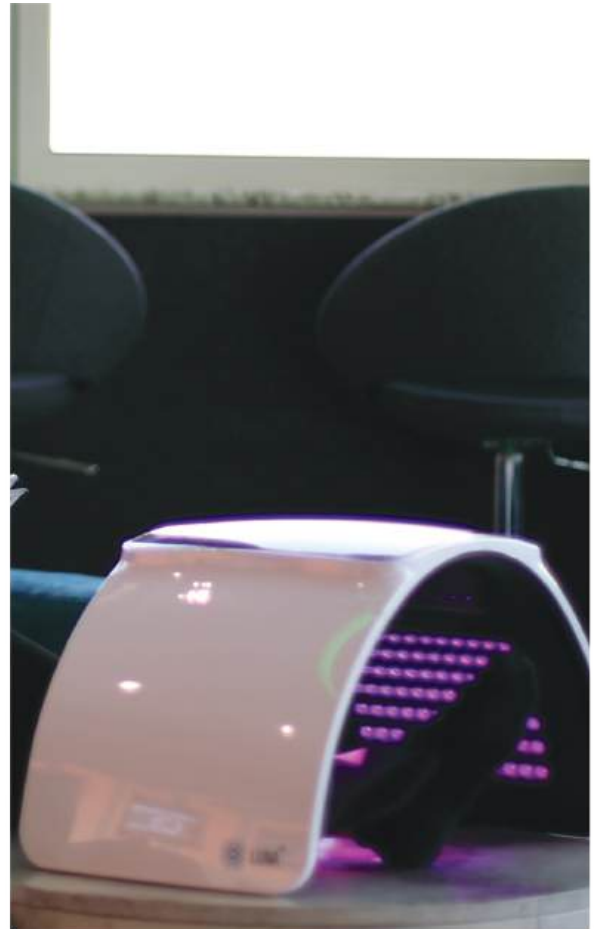

USING THE DEVICE

To select your treatment, start by setting the duration of your treatment by adjusting the plus and minus icons. You can increase or decrease the timer in increments of 5 minutes.

NEXT, SELECT YOUR LED LIGHT THERAPY MODE:

RED LIGHT TREATMENT:

To initiate the red light mode, tap the icon named, "RED". To deselect the red mode, tap the icon again.

BLUE LIGHT TREATMENT:

To initiate the blue light mode, tap the icon named, "BLUE". To deselect the blue mode, tap the icon again.

GREEN LIGHT TREATMENT:

To initiate the green light mode, tap the icon named, "GREEN". To deselect green mode, tap the icon again.

CYCLE TREATMENT (ALL LIGHTS COMBINED):

To initiate the cycle mode, tap the icon named, "CYCLE". To deselect the cycle mode, tap the icon again.

Other treatments (orange, yellow, and violet light):

To initiate any other light modes including orange, yellow, or purple, tap the icon named, "OTHERS". Press the icon once to initiate the orange mode, twice to initiate the yellow mode, and thrice to initiate the purple mode. To deselect the mode, press the corresponding button again, the same number of times.

CALCIUM SUPPLEMENTATION IRRADIATION MODE:

is an extra treatment that can be added to any of the other available treatments or done by itself. Tap the CA2+ icon to enter this mode. The default time set for this mode is 30 minutes. However, you can increase or decrease the duration of this treatment by adjusting the plus and minus icons. To deselect the calcium supplementation mode, tap the CA2+ icon again. Unless prescribed by your medical provider, we highly recommend that you do not exceed 30 minutes of the calcium supplementation irradiation treatment per day.

LOCK/UNLOCK

Once you have selected your desired treatment, tap the lock icon with a single short press to set the treatment. To unlock the treatment press and hold the lock icon for 3 seconds. You are now ready to place the device above the treatment area and start your treatment! End each treatment with your favorite moisturizer or serum for maximum results.

IMPORTANT

When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.





SPACETOUCH

TREATMENT INDEX

Psoriasis	16	Eczema	24
Neck Pain	17	Dull Skin	25
Hair Growth	18	Sun Damage & Solar Keratoses (sunspots)	26
Cellulitis	19	Pigmentation	27
Sagging Skin	20	Acne	28
Rosacea	21	Wrinkles	29
Scars	22	Asthma	30
Stretch Marks & Spider Veins	23	Low Bone Density & Osteoporosis	31

PSORIASIS

ROLANDE FRANKE



TARGET
Psoriasis



AREAS OF USE
Face, body, and all
areas of concern



DOSE
35 Minutes
per treatment



TREATMENT LENGTH
3 Times per week
for 3 months

DIRECTION

This treatment utilizes red LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 35 minutes, and select the red mode. Place the device above the desired treatment area for 35 minutes. Once the red LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.

 **TIMER**
35 Minutes

 **LIGHT**
Red



NECK PAIN

SEAN MERLING



TARGET

Neck pain



AREAS OF USE

Neck



DOSE

60 Minutes
per treatment



TREATMENT LENGTH

3 Times per week
for 3 months

DIRECTION

This treatment utilizes red LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 60 minutes, and select the red mode. Lie flat on your stomach, and place the device around your neck for 60 minutes. Once the red LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER

60 Minutes



LIGHT

Red



HAIR GROWTH

ROLANDE FRANKE



TARGET

Hair growth



AREAS OF USE

Head



DOSE

30 Minutes
per treatment



TREATMENT LENGTH

3 Times per week
for 3 months

DIRECTION

This treatment utilizes red LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 30 minutes, and select the red mode. Lie down and place the device around your scalp for 30 minutes. Once the red LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER

30 Minutes



LIGHT

Red



CELLULITIS

DR. JANA BELVIEW



TARGET
Cellulitis



AREAS OF USE
Face, body, and all
areas of concern



DOSE
30 Minutes
per treatment



TREATMENT LENGTH
5 Times per week
for 3 months

DIRECTION

This treatment utilizes red LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 30 minutes, and select the red mode. Place the device above the desired treatment area for 30 minutes. Once the red LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER
30 Minutes



LIGHT
Red



SAGGING SKIN

DR. JANA BELVIEW



TARGET

Sagging skin



AREAS OF USE

Face, body, and all areas of concern



DOSE

75 Minutes
per treatment



TREATMENT LENGTH

2 Times per week
for 3 months

DIRECTION

This treatment utilizes green LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 75 minutes, and select the green mode. Place the device above the desired treatment area for 75 minutes. Once the green LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER

75 Minutes



LIGHT

Green



ROSACEA

DR. MICHAEL MCGEE



TARGET
Rosacea



AREAS OF USE
Face, body, and all
areas of concern



DOSE
45 Minutes
per treatment



TREATMENT LENGTH
4 Times per week
for 3 months

DIRECTION

This treatment utilizes yellow LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 45 minutes, and select the yellow mode. Place the device above the desired treatment area for 45 minutes. Once the yellow LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.

 **TIMER**
45 Minutes

 **LIGHT**
Yellow



SCARS

SEAN MERLING



TARGET
Scars



AREAS OF USE
Face, body, and all
areas of concern



DOSE
60 Minutes
per treatment



TREATMENT LENGTH
4 Times per week
for 3 months

DIRECTION

This treatment utilizes violet LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 60 minutes, and select the violet mode. Place the device above the desired treatment area for 60 minutes. Once the violet LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER
60 Minutes



LIGHT
Violet



STRETCH MARKS & SPIDER VEINS

DR. MICHAEL MCGEE



TARGET

Stretch marks & spider veins



AREAS OF USE

Body and all areas of concern



DOSE

60 Minutes per treatment



TREATMENT LENGTH

4 Times per week for 3 months

DIRECTION

This treatment utilizes orange LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 60 minutes, and select the orange mode. Place the device above the desired treatment area for 60 minutes. Once the orange LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER

60 Minutes



LIGHT

Orange



ECZEMA

DR. MICHAEL MCGEE



TARGET
Eczema



AREAS OF USE
Face, body, and all
areas of concern



DOSE
40 Minutes
per treatment



TREATMENT LENGTH
3 Times per week
for 3 months

DIRECTION

This treatment utilizes orange LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 40 minutes, and select the orange mode. Place the device above the desired treatment area for 40 minutes. Once the orange LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER
40 Minutes



LIGHT
Orange



DULL SKIN

DR. MICHAEL MCGEE



TARGET
Dull skin



AREAS OF USE
Face, body, and all
areas of concern



DOSE
40 Minutes
per treatment



TREATMENT LENGTH
3 Times per week
for 3 months

DIRECTION

This treatment utilizes orange LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 40 minutes, and select the orange mode. Place the device above the desired treatment area for 40 minutes. Once the orange LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.

 **TIMER**
40 Minutes

 **LIGHT**
Orange



SUN DAMAGE & SOLAR KERATOSIS (SUNSPOTS)

ROLANDE FRANKE



TARGET

Sun damage and sunspots



AREAS OF USE

Face, body, and all areas of concern



DOSE

60 Minutes per treatment



TREATMENT LENGTH

5 Times per week for 3 months

DIRECTION

This treatment contains two steps, one using the yellow LED light and the other using the blue LED light. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 40 minutes, and select the yellow mode. Place the device above the desired treatment area for 40 minutes. Once the yellow LED light treatment is complete, set the timer for 20 minutes and select the blue mode. Place the device above the same treatment area for 20 minutes. After completing all treatments, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER

40 & 20 minutes



LIGHT

Yellow & blue



PIGMENTATION

DR. JANA BELVIEW



TARGET

Pigmentation



AREAS OF USE

Face, body, and all areas of concern



DOSE

60 Minutes per treatment



TREATMENT LENGTH

3 Times per week for 3 months

DIRECTION

This treatment utilizes both blue and green LED lights. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 60 minutes, and select the blue and green modes at the same time. Place the device above the desired treatment area for 60 minutes. Once treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER

60 Minutes



LIGHT

Blue and green



ACNE

ROLANDE FRANK



TARGET
Acne



AREAS OF USE
Face, body, and all
areas of concern



DOSE
45 Minutes
per treatment



TREATMENT LENGTH
4 Times per week
for 3 months

DIRECTION

This treatment utilizes both red and blue LED lights. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 45 minutes, and select the red and blue modes at the same time. Place the device above the desired treatment area for 45 minutes. Once treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER
45 Minutes



LIGHT
Red & Blue



WRINKLES

DR. JANA BELVIEW



TARGET
Wrinkles



AREAS OF USE
Face & neck



DOSE
60 Minutes
per treatment



TREATMENT LENGTH
3 Times per week
for 6 months

DIRECTION

This treatment utilizes both red and blue LED lights. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 60 minutes, and select the red and blue modes at the same time. Place the device above the desired treatment area for 60 minutes. Once treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.

 **TIMER**
60 Minutes

 **LIGHT**
Red & blue



ASTHMA

DR. ACHIM LEUTH



TARGET
Asthma



AREAS OF USE
Chest and neck



DOSE
20 Minutes
per treatment



TREATMENT LENGTH
3 Times per week
for 4 months

DIRECTION

This treatment utilizes green LED light therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 20 minutes, and select the green mode. Place the device above the desired treatment area for 20 minutes. Once the green LED light treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.

 **TIMER**
20 Minutes

 **LIGHT**
Green



LOW BONE DENSITY & OSTEOPOROSIS

DR. ACHIM LEUTH



TARGET

Low bone density & osteoporosis



AREAS OF USE

Face & neck



DOSE

30 Minutes
per treatment



TREATMENT LENGTH

30 Times per week
for 6 months

DIRECTION

This treatment utilizes Calcium supplementation irradiation therapy. Connect Luna to a power outlet and place it on a safe, flat surface. Tap the power icon to turn on the device. Set the timer to 30 minutes, and select the CA2+ mode. Place the device above the desired treatment area for 30 minutes. Once treatment is complete, apply your favorite cream or one recommended by a specialist.

IMPORTANT: When placing the device over the facial area, it is essential to wear the protective glasses that come with the device to avoid damage to your eyes.



TIMER

30 Minutes



TREATMENT

Calcium supplementation irradiation



CONTRAINDICATIONS, WARNINGS, AND PRECAUTIONS

CONTRAINDICATIONS

DO NOT use the device if you are pregnant

DO NOT use the device while showering, bathing or swimming

DO NOT drop the device or subject it to excessive force

DO NOT continue to use the device if you develop a skin irritation; consult a doctor if skin irritation persists

DO NOT use the device if you suffer from lupus erythematosus, albinism, photosensitive disorders, or any medications that can cause photosensitivity unless you have consulted with a doctor. LED lights can cause severe skin reactions for the conditions listed

DO NOT use the device if you suspect that it is damaged or needs repair. Instead, contact Spacetouch for a replacement

WARNINGS

DO NOT use the device without removing all makeup prior to treatment. Some ingredients in these products may cause light sensitization or reduce the efficiency of the treatment

DO NOT tamper with or attempt to modify the device in any way. Modifying any Spacetouch device will result in a voided warranty

DO NOT continue using the device if it begins to overheat or feel uncomfortably warm on your skin, turn off it immediately and contact customer support

DO NOT continue using the device if you encounter prolonged skin redness or discomfort that lasts for more than 2 hours after treatment. This could be due to light sensitivity or another condition. Please consult your doctor for advice

PRECAUTIONS

Before using Luna, we highly recommend that you conduct a light sensitization test. Stop using the device if you experience any adverse effects from the treatment and consult with your doctor if you have any concerns about the product and/or treatment.

LIGHT SENSITIVITY TEST

Step 1: Power the device and use it on a clean face.

Step 2: Use the device for one round of treatment. The device will turn off once the treatment is complete. At any point, if you experience uncomfortable heat or discomfort, stop using the device immediately.

Step 3: Wait for 24 hours before using the device again. If you experience pain, discomfort, or redness that lasts for more than 2 hours, please consult a doctor for advice as you might have light sensitivity.

LIMITED LIFETIME WARRANTY

The lifetime warranty for Spacetouch devices guarantees coverage for manufacturing defects. To activate the warranty, you must register your device on our website, www.spacetouch.com, within 7 days of receiving the order. Failure to register the device correctly and timely will invalidate the warranty. Please note that the warranty does not cover damages caused by intentional or unintended misuse, use of the device for purposes other than its intended use, or any damage beyond normal wear and tear. Certain restrictions may apply.

REGISTER YOUR DEVICE BY
SCANNING THE QR CODE BELOW



SPACETOUCH



@spacetouch

