

SPACETOUCH

ADVANCED SKINCARE THERAPY

REVIVER



spacetouch.com

REVIVER USER GUIDE

Hello!

Thank you for purchasing Reviver and welcome to the Spacetouch family! There is no better investment than investing in yourself! You are now part of a growing community of millions of people who use advanced LED skincare therapy to achieve flawless skin.

As you unbox your new device, take a few moments to review the manual and become familiar with its features. The instructions provide important details on how to safely operate your device.


Spacetouch is committed to providing you with the best skincare experience possible. If you have any questions or concerns, our friendly and knowledgeable customer service team is here to help. You can reach us at Service@Spacetouch.com.


TABLE OF CONTENTS

Configuration	03
Welcome	05
The Heads	06
Ampoule Series	10
Adding or Removing Ampoules	12
Operating The Device	15
Maintenance and Storage Tips	17
Charging Your Device	17
Cleaning and Maintenance Guidelines	19
User Safety Guidelines	20
Device Safety Guidelines	23
Product Specifications	26
Frequently Asked Questions	27
Safe Disposal and Care Guidelines	28
Device Registration	30


CONFIGURATION

Body	X1
IR+Red & Topical Heat Head	X1
EMS Head	X1
Microneedle Head	X1
3ml Ampoule Bottle	X12
Charging Base	X1
Needle	X1
Charging Cable	X1

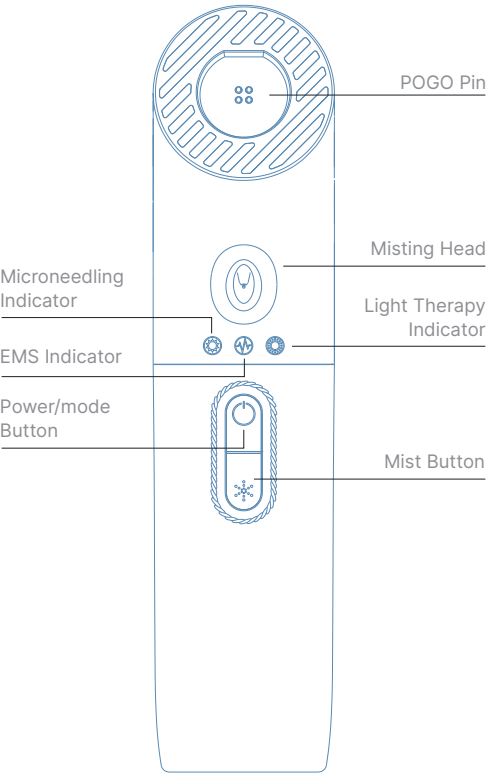

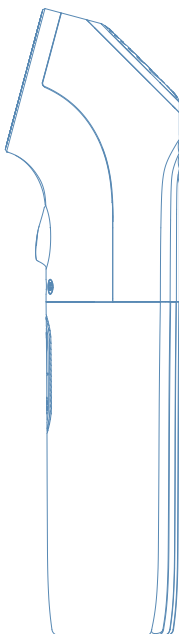


Side

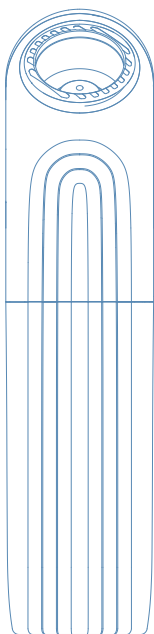

back



EMS Head

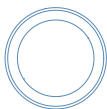

Light Therapy Head



Microneedling Head


Ampoule Bottle 3ml



Silicone Ring



Charging Base

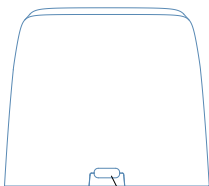

Type-C
Charging Port

Welcome To The Future of Skincare With Your New REVIVER Device.

Embark on a revolutionary journey towards flawless skin with this state-of-the-art device, meticulously engineered to transform your skincare routine. Reviver isn't just a device; it's a compact powerhouse, blending sleek design with cutting-edge functionality. It comes equipped with three specialized treatment heads, each tailored for different skincare needs, ensuring a comprehensive approach to skin health.

What truly sets Reviver apart is its innovative nano serum mister. This feature, incorporating patented technology, is designed to deeply penetrate your skin with high-grade serum concentrates, enhancing the efficacy of each treatment.

Our mission in designing Reviver was to bring the sophistication of professional skincare tools into your hands. Whether you're a seasoned professional or a home skincare enthusiast, this device breaks down the barriers, making advanced skincare accessible and convenient for everyone. It's crafted for simplicity, efficiency, and portability so that you can enjoy a professional skincare experience anywhere, anytime.

Four years of rigorous research and development have culminated in Reviver,


our most compact yet advanced device to date, with seven patents pending. Its robust design includes easily interchangeable magnetic heads, allowing you to customize your skincare regimen. The three included magnetic treatment heads are the foundation of your journey towards professional-level skincare right in the comfort of your home.

Embrace the elegance and effectiveness of Reviver, and step into a world where premium skincare is at your fingertips.

THE HEADS

SKIN REPAIR

Microneedling Head

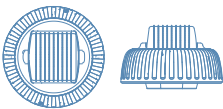

Dive into the world of microneedling, a groundbreaking technique designed to rejuvenate your skin. This method is renowned for its ability to significantly enhance your skin's appearance and texture. Whether you're looking to reduce stretch marks, diminish scars, smooth out wrinkles, or address other skin imperfections, microneedling is your go-to solution.

Unlike certain facial rejuvenation treatments that rely on heat, or lasers,


microneedling stands out with its minimal risk of skin discoloration. It's a safe and effective way to achieve a more youthful, radiant complexion.

Usage Recommendation:

- For best results, we recommend using the microneedling head **once every 4 to 6 weeks on clean skin for 15 minutes each session.**
- Pair the microneedling process with one of our specially formulated nano mist ampoules for an even more remarkable outcome. This combination is designed to amplify the benefits and leave your skin feeling revitalized and glowing.

ANTI-AGING, HYDRATING, AND COLLAGEN PRODUCTION

Infrared (IR)+RED & Topical Head

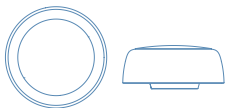

Discover the secret to youthful, radiant skin with our cutting-edge IR+RED Light Therapy Head. Harnessing the power of both infrared and red light, this innovative tool is a favorite among skincare experts and medical professionals. It's designed to visibly diminish wrinkles and fine lines, soothe inflammation, accelerate healing

post-procedures, and rejuvenate your skin.

What makes this therapy so effective? It's all about enhancing your skin's cellular health. Red light therapy boosts collagen production, a key ingredient for firm, smooth skin. You'll notice improvements in your skin's tone, texture, and elasticity.

Usage Recommendation:

- For best results, we recommend using the IR+RED and Topical Heat Head **2 to 5 times weekly on clean skin** for 15 minutes each session.
- For amplified benefits, pair each session with one of our exclusive nano mist ampoules.

LIFTING, TONING, & BRIGHTENING

EMS Head



Elevate your skincare routine with our advanced EMS (Electrical Muscle Stimulation) technology head. Designed to rejuvenate your skin from within, EMS works wonders by gently lifting and toning your skin. It boosts the production of collagen and elastin, essential elements your skin needs to rebuild its inner structure. As a result,

you'll notice a remarkable increase in the elasticity and thickness of your skin.

Not only does EMS enhance your skin's texture and firmness, but it also significantly diminishes the appearance of cellulite, sagging skin, and stretch marks. It's a revolutionary way to achieve firmer skin without surgery or injections.

This head also features cutting-edge Blue LED light therapy. Utilizing specific wavelengths of blue light (405-410 nanometers), this therapy delves deep into your skin to eradicate bacteria. It's exceptionally effective in reducing skin issues like pigmentation and discoloration, leaving your skin clearer and more radiant.

By combining the prowess of EMS and Blue LED therapy, this head is a powerhouse for facial toning and lifting. It works synergistically to brighten your skin and bestow it with a flawless glow.

Usage Recommendation:

- For best results, we recommend using the EMS and Blue LED Light Head **2 to 5 times weekly on clean skin for 15 minutes each session.**
- Enhance your experience by pairing it with one of our exclusive nano mist ampoules, designed to amplify the

benefits and leave your skin looking it's absolute best.

THE SPACETOUCH AMPOULE SERIES: TAILORED FOR FACE & BODY

Our Spacetouch Ampoule range is a revolutionary approach to skincare, offering both facial and body treatments.

Utilizing cutting-edge nano mist technology, the ampoules transform into a fine mist, ensuring deep penetration and optimal absorption into your skin's pores. This process amplifies the effectiveness of the active ingredients, providing remarkable results.

THE AMPOULES

Experience the power of precision skincare with our specially formulated ampoules.

Each tiny vial is packed with potent, high-concentration ingredients, crafted to target specific skin concerns with pinpoint accuracy.

Whether you're looking to hydrate, brighten, or rejuvenate, these ampoules deliver just the right dose, right where you need it.

EXPLORE OUR SIGNATURE FORMULAS

Your Reviver device comes with an

exclusive selection of four sample formulas, each designed to address a different skin need:

- Repairing - For skin that needs a healing touch
- Hydrating & Glowing - To infuse your skin with moisture and a radiant sheen
- Lifting & Toning - For a firmer, more sculpted appearance
- Brightening - To illuminate and even out your skin tone

EFFORTLESS APPLICATION WITH THE REVIVER DEVICE

The Reviver device puts control in your hands. Adjust the amount of ampoule used by simply pressing the mist button. This user-friendly feature allows you to evenly distribute the mist over your skin, ensuring a uniform application every time.

STAY STOCKED UP

To replenish your ampoules or explore other formulas, visit our website at spacetouch.com.

Alternatively, you can consult a skin specialist at your nearest skincare clinic for personalized recommendations. Let your skin experience the future of

skincare with Spacetouch.


ADDING OR REMOVING AMPOULES

INSERTING AN AMPOULE


- 1 Prepare the Ampoule:** Each ampoule contains a precise 3ml dose. Begin by carefully removing the cap to open the ampoule. Ensure the ampoule is kept upright during this process.



- 2 Attach the Ampoule:** Gently slide the ampoule into the device, aligning it with the Silicon protector. Remember, the ampoule should remain upright – do not invert it. Instead, invert the device body, guiding the ampoule into position until it clicks securely into place.

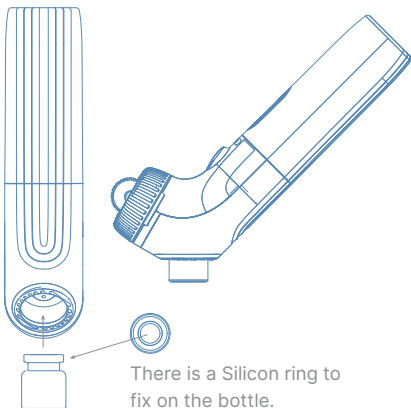


- 3 Ready to Use:** Once the ampoule is firmly attached, you can return the device body to its upright position. It's now ready for use.

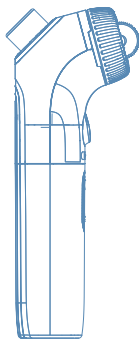


REMOVING OR REPLACING AN AMPOULE

INSERTING AN AMPOULE


- 1 Steady Positioning:** Hold the device so that the ampoule is upright, just as you did when inserting it.

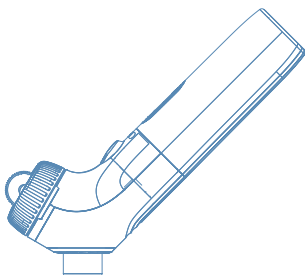

- 2 Gentle Removal:** Firmly grasp the ampoule and carefully lift the device upwards, detaching the ampoule. If there's any remaining serum, ensure it drains back into the ampoule.

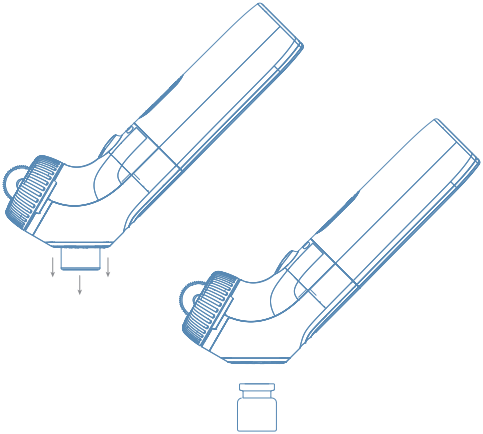

- 3 Clean Up:** If necessary, wipe the ampoule feeder with a dry cloth to maintain cleanliness.

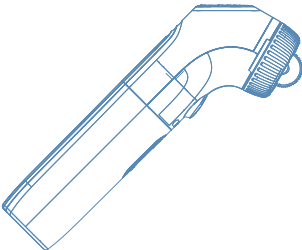

OPERATING THE DEVICE

Your Reviver device features a state-of-the-art automatic recognition head system.

This innovative design ensures that all buttons on the device can be used interchangeably without any operational issues.

The machine is programmed to automatically shut down after 30 minutes of continuous operation for your convenience and safety.

USING THE MISTER

To Activate the Mister:

- 1** Press the Mist button.
- 2** To deactivate, press the Mist button.

MICRONEEDLE HEAD

Setting Up the Microneedle Head:

- 1** Attach the Microneedle Head onto the device. The Microneedle indicator light will turn on, signaling readiness.
- 2** Apply the serum with the sprayer before or while using the Micro Needle roller head for optimal results.

EMS (ELECTRICAL MUSCLE STIMULATION) HEAD

EMS Added 3 Gears

Activating the EMS Mode:

- 1** Attach the EMS head to the device. It will automatically switch to EMS mode. Short Press the power button again go adjust the gear, there are 3 gears.
- 2** Press the "Power" button for about 1.5 seconds to activate the EMS and blue light mode. The EMS indicator light will signal that the mode is active.
- 3** Press and hold the power button for about 1.5 seconds again to turn off the EMS mode. The EMS indicator light will turn off.
- 4** For amplified benefits, pair each session with one of our exclusive nano mist ampoules.

IR+RED & TOPICAL HEAT HEAD

Engaging Red Light Mode:

- 1** Attach the Red Light Head onto the device. It will automatically enter red light mode.
- 2** Press the Power button for approximately 1.5 seconds to activate the red light therapy. The light therapy indicator light will illuminate to show it's active.

Press and hold the Power button for

around 1.5 seconds again to deactivate the therapy mode. The light therapy indicator light will switch off.

MAINTENANCE AND STORAGE TIPS

- 1** Always place your device back on its charging base after each use. This ensures it's always charged and ready for your next session.
- 2** Store your ampoules in a cool, dry place. A room or cabinet away from moisture works best. Remember, there's no need to refrigerate the ampoules.
- 3** Keep the heads of your device clean to ensure optimal performance. Regular cleaning also extends the life of your device.
- 4** The device body loves to stay dry! It's not waterproof, so remember not to take it into the shower or wash it under running water. A damp cloth for the exterior is all you need for cleaning.

CHARGING YOUR DEVICE

- **Low Power Alert:** If your device is about to run out of power, the power indicator will blink 5 times before the device automatically powers down.

- **Required Adapter Specifications:**
Use an adapter with an output voltage of 5V DC.
- **Charging Options:** You have two ways to charge your device:
 - Direct Cable Charging: Connect your device to a power source using the provided cable.
 - Base Charging: Place your device on the charging base. When the base is connected to power, its indicator light will turn on.
- **Charging Indicators:**
 - During Charging: While your device is charging, the power light on the device will flash, signaling that charging is in progress.
 - Full Charge: It typically takes 2 to 3 hours to fully charge the device. You'll know it's fully charged when the power light turns off.
- **Usage During Charging:** To ensure safety and device integrity, your device cannot be used while charging.

CLEANING AND MAINTENANCE GUIDELINES

DEVICE BODY

- The device body is not waterproof.
- AVOID spraying or immersing it in water.
- To clean the surface, gently wipe it down with a damp cloth or wet tissue.

MISTER

- Clean the sprayer using either distilled water or medical alcohol for optimal performance.
- After refilling, spray out any residual water.

MICRONEEDLE HEAD ROLLER HEAD

- The Microneedle head is waterproof and can be cleaned or disinfected by immersing it in medical alcohol or water.
- Handle the micro needle with care to avoid damage.
- After each use, clean and dry the

Microneedle head thoroughly.
Store it in a dry place, away from
direct light to maintain its condition.

EMS HEAD

- Clean the head using either distilled water or medical alcohol for optimal performance.
- Do not use harsh chemicals or corrosive liquids for cleaning, as they can damage the head.

RED LIGHT HEAD

- Clean the Red light head using water.
- Similar to the EMS head, avoid using corrosive liquids to clean the Red light head to preserve its functionality.

USER SAFETY GUIDELINES

- 1 Important:** Do not use the device near medical devices like pacemakers, heart-lung machines, electrocardiographs, or IV drips.
- 2 Gentle Usage:** Avoid lingering on one area of the skin for too long or overusing the device on the same spot.
- 3 Post-Use Skin Reaction:** It's normal to experience slight skin irritation or redness after use, as the product

stimulates skin circulation. These reactions are temporary and should fade on their own.

- 4 Skin Care Precautions:** Do not use this product on areas recently treated with cosmetic procedures such as plastic surgery, rhinoplasty, or hair removal. Avoid using it if you're feeling unwell, overly fatigued, or if you notice unusual skin conditions. It's also not suitable for use on the eyes, eyelids, or lips.
- 5 Consult Healthcare Providers:** This is a beauty care product, not intended for medical use. If you're under medical treatment, please consult your doctor before use, especially if you have a serious illness, have recently received treatment requiring rest, or suffer from acute illnesses.
- 6 Initial Usage:** You may notice small wavelet dots in a new device; these are simply signs of our thorough quality check. Feel confident using it!
- 7 Handling and Storage:** Please do not disassemble, repair, or modify this product. Keep it away from flammable materials and out of childrens reach. It's not designed to be water-resistant, so avoid using it in damp environments like bathrooms or while showering.

- 8 Cleaning Instructions:** After use, gently wipe the device body with a damp cloth (avoid using harsh solvents like gasoline, thinner, or alcohol). Be sure your hands are dry when plugging or unplugging the adapter. Regularly clean the power cord to prevent damage and dust accumulation.
- 9 Careful Handling:** Avoid dropping or subjecting the product to impact. Before use, check for any signs of damage, such as cracks or deformities. Stop using it immediately if you notice any damage or malfunction.
- 10 Remove Jewelry:** Remove any face or neck accessories like earrings to ensure a safe and comfortable experience.
- 11 Responsibility Clause:** Failure to adhere to these guidelines absolves us of any responsibility for injuries, malfunctions, or damage resulting from incorrect usage, disassembly, or unauthorized modifications.

DEVICE SAFETY GUIDELINES

MISTER GUIDELINES

- **Particle-Free Liquids:** Only use original Spacetouch Ampoules to prevent clogging.
- **Cleaning:** Rinse the sprayer with clean water before and after each use.
- **Unclogging:** If the nozzle becomes blocked, gently unblock it using the provided needle.
- **Skin Safety:** Avoid using the sprayer on damaged, ulcerated, or inflamed skin.

MICRONEEDLE HEAD GUIDELINES

- **Preparation:** Always clean the device before using the Microneedle Head.
- **Disinfection:** Sterilize the head with medical alcohol before use to prevent bacterial infections.
- **Skin Preparation:** Ensure your skin is clean, removing any dirt or oil.
- **Application Technique:** Serum can be applied before or during rolling.

- **Gentle Application:** Roll gently without excessive pressure to avoid skin damage.
- **Sensitive Areas:** Avoid sensitive skin areas like eyes and lips.
- **Post-Use Care:** Apply a medical facial mask after rolling to aid in skin recovery.
- **Maintenance:** Clean and dry the roller head after each use. Store it in a dry place away from direct light.
- **Roller Head Replacement:** Replace the Microneedle head every Five to Six months.
- **Safety Note:** Do not use the roller if the needles are deformed, contaminated, or if the skin is damaged, ulcerated, or inflamed.

EMS HEAD GUIDELINES

- **Combination Use:** The EMS head can be used in conjunction with the red light feature.
- **Sequence of Use:** Utilize the EMS head before, but not after, Microneedling .
- **Conductive Medium:** Always use the

EMS head with a Spacetouch Ampoules.

- **Skin Safety:** Avoid using the EMS head on damaged, ulcerated, or inflamed skin.

RED LIGHT HEAD

- **Combination Use:** The red light head can be used alongside the EMS head.
- **Flexible Application:** Feel free to use the red light feature either before or after Microneedling.
- **Temperature Caution:** The red light head operates at a standard temperature of 42°C (107.6°F). Discontinue use immediately if you experience any burning sensation.
- **Skin Safety:** Do not use the red light head on skin that is damaged, ulcerated, or inflamed.

PRODUCT SPECIFICATIONS



PARAMETERS

Product size	75*69*198mm
Rated power	12W
Battery capacity	1200mAh
Voltage	5V
Charging time	2~3hours
Needle length	0.5mm
Needle material	SUS304
Red light wavelength	630nm
Infrared light wavelength	840nm
Tank capacity	3ml
Rated Temperature	42 °C

FREQUENTLY ASKED QUESTIONS

1. The device does not automatically recognize the corresponding functional head.

What To Do: Ensure the conductive spring pin of the device is popping out as it should. Also, check if the contact points on the functional head are clean and free from dirt.

2. Is it okay to use thick serum?

Do not use any other serums aside from Spacetouch Ampoules.

3. Does the Microneedle head cause skin pain during use?

A slight pricking sensation is normal when using the Microneedle Head. If it doesn't roll smoothly, stop immediately to prevent skin scratches. Remember, never drag the roller sideways.

5. The device won't turn on, or power indicator is blinking?

This usually indicates a low battery. Connect the device to a DC 5V power source for charging. Typical charging time ranges from 3 to 4 hours.

SAFE DISPOSAL AND CARE GUIDELINES

RECYCLE AND CARE FOR YOUR RECHARGEABLE BATTERIES (LITHIUM-ION)

To maintain safety and support environmental efforts, please follow these guidelines for your device rechargeable battery:

- 1 Battery Exclusivity:** Use the device's battery only for this product. Avoid charging a battery that's been removed.
- 2 Handling the Device body:** Don't disassemble, hit, squeeze, or expose the device to fire.
- 3 Swelling Indicator:** Stop using the device if noticeable swelling occurs.
- 4 Temperature Precautions:** Avoid exposing the device to high temperatures.
- 5 Wet Battery Protocol:** Stop usage if the battery becomes wet.
- 6 Physical Damage:** Do not puncture, forcefully hit, or alter the device's structure.
- 7 Battery Safety:** Avoid directly

connecting the battery's positive and negative terminals with metal.

- 8 **Storage Precautions:** Keep it away from small metal objects like necklaces or hairpins.
- 9 **Child Safety:** Keep the battery away from infants to prevent accidental ingestion. In case of ingestion, seek medical help immediately.

10 Leak Response:

- a **Eye Contact:** If battery fluid contacts your eyes, avoid rubbing them. Rinse with water and seek immediate medical help to prevent potential blindness.
- b **Skin Contact:** If battery fluid touches your skin or clothes, it may cause irritation or injury. Rinse affected areas with water and get medical assistance if discomfort persists.

- 8 **Disposal and Repair:** The device should only be disassembled when being disposed of to prevent risks like fire or electric shock. For repairs (excluding battery replacements), please contact the store of purchase.

Battery Disposal: When disposing of the device, first turn it off and remove the lithium-ion battery. Hand over the battery to a professional

recycling agency specializing in
rechargeable batteries.

Register your device by
scanning the QR code below





@spacetouch